



Small Animal Decapitator


Cat. No. 7950

General

This Decapitator is a simple, rugged, hand operated shear for humane, instant decapitation of laboratory animals.

With its 5 cm square maximum blade opening, it is suitable for decapitating guinea pigs, rats, mice and other small animals.

The blades are made of stainless steel, ground and induction hardened to Rockwell 60, to re-tain a razor sharp edge.

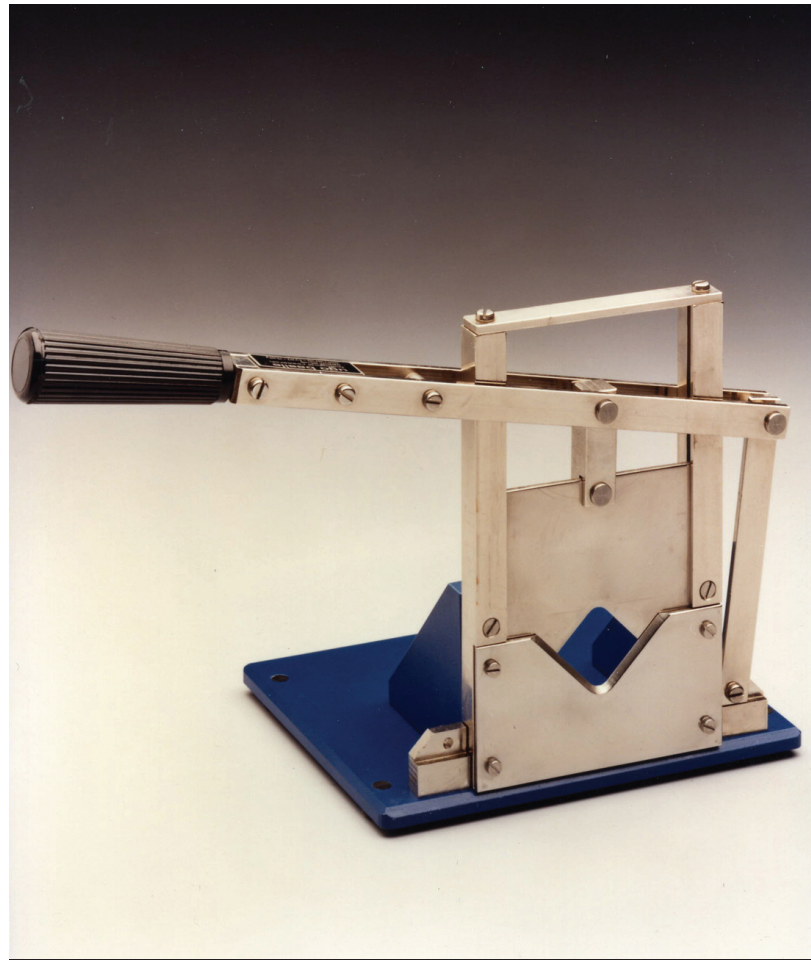

for humane, instant decapitation
of small laboratory animals

Main Features

- 5 cm square opening
- Induction-hardened stainless steel blades
- Heavy gauge aluminium base with four mounting holes

General

The 7950 decapitator cuts cleanly and quickly through bone and tissue.

Replacement blades are available and must be installed at the factory.

The heavy gauge aluminum base has four holes for bench mounting.

The base is coated with textured oven-cured ep-oxy paint; stud bolts, screws and washers are all made of stainless steel; the frame, the lever handle and the linkage of the upper blade are made of nickel plated brass.

A unique feature allows the handle to be re-located for right or left handed use. Please specify your preference when ordering, by adding -R (right) or -L (left) to the catalogue number.

The unit can be totally dipped for cleaning.

Ordering Information

7950 Small Animal Decapitator

7951 Spare Blade set, two pieces

NOTE blades must be installed/replaced at our factory.

