



Ambulatory EEG recordings are undoubtedly the best way to:

- Determine the presence of pseudo seizures
- Determine the presence of Absence seizures
- Aid in titrating anticonvulsants
- Diagnose more patients in your practice at the lowest cost

Ambulatory EEG has been an underutilized diagnostic tool because:

- High cost of equipment
- Poor quality of recordings
- Available products had a limited number of channels

The new Compumedics **Safiro** system pushes aside many of those historical limitations. **Safiro** offers the following advantages, many never before available for AEEG recordings:

Convenience

- Only 9.6oz (300gms) with internal rechargeable batteries, the lightest recorder on the market
- Only 139.7mm x 76.2mm x 38.1mm, 5.5"x 3"x 1.5", the smallest recorder on the market
- Batteries recharge in the recorder; no need to open Safiro between recordings. Just plug it in.
- Records on tiny compact flash card, up to 1GB

Capability

- Record up to 32 referential channels
- Record at sampling rates up to 256Hz
- Record at resolutions up to 16 bit
- Safiro is capable of multi-day ambulatory recordings

Confidence

- Wireless, infrared (IR) link for pre-check impedance measurements and signal calibration, data stored with the patient recording
- See all channels simultaneously during set up on any notebook or desktop computer, just as you do in the laboratory

safIRO™


Ambulatory EEG



Set up your patient and select the montage - up to 32 channels.



Observe the EEG data, calibration signals and impedance values. Save the impedance values to be reviewed with the final recording.

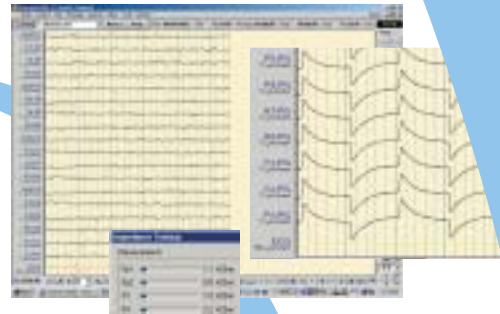

Compumedics' superior digital technology allows simple operation and better quality recordings



Send the patient home.

Review the data and produce the report with ProFusion EEG software. Spike and Seizure detection software is optional.



When the recorder is returned the next day, just insert the flash card into the computer. The recording is downloaded quickly and automatically.

safIROTM
Ambulatory EEG

ProFusionTM
EEG
EEG Analysis Software

For more information please contact:

Compumedics USA, Inc:
7850 Paseo del Norte
Suite 101
El Paso, Texas 79912
Toll Free: 1 877 717 3975
Fax: 1 915 845 2965

Compumedics Limited, Australia:
30-40 Flockhart Street
Abbotsford VIC 3067
Australia
Tel: +61 3 8420 7300
Fax: +61 3 8420 7399
Free Call: 1800 651 751

Compumedics Singapore, Pte Ltd:
7500A Beach Road,
#10-323 The Plaza,
Singapore 199591
Tel: +65 6392 0902
Fax: +65 6392 0570

Compumedics Germany, GmbH:
Heussweg 25
20255 Hamburg Germany
Tel: +49 40 40 18 99 41
Fax: +49 40 40 18 99 49



www.compumedics.com

