



Anesthetizing Box

Cat. No. 7900 (rodents) 7910 (rabbits)

General

Our Induction Boxes are conveniently dimensioned induction boxes, featuring a sliding lid. They are made of Perspex and prove to be particularly useful to confine laboratory animals during anesthetizing.

The transparent acrylics permits the animal to be kept under constant observation.


Two tubing connectors of nickel plated brass are fitted into each end, one located at the top of the box and the other at the bottom.

Any (non-explosive!) gas mixture can be used. In case small quantities of liquid, as ether or chloroform are used, soak a cotton wool flock and place it in a small Becker, in-side the box.



Our Induction chambers are ideal to work with our new Anesthesia Systems

TO CONFINE SMALL LABORATORY ANIMALS DURING ANESTHETIZING

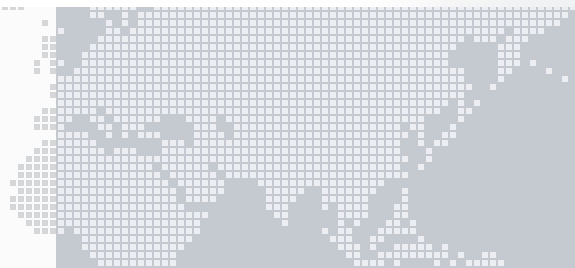

Box Dimensions:

- | | |
|------------------------------------------|--------------|
| ● 7900 - induction box for small rodents | 25x13x13 (h) |
| ● 7910 - induction box for rabbits | 40x22x21 (h) |